
Munddusche

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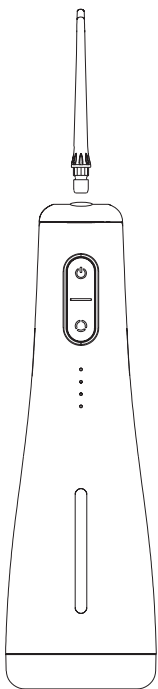
Water flosser

INSTRUCTION MANUAL

Read these instructions carefully before using the device and keep it safe for future reference.

KLAMER

Water flosser

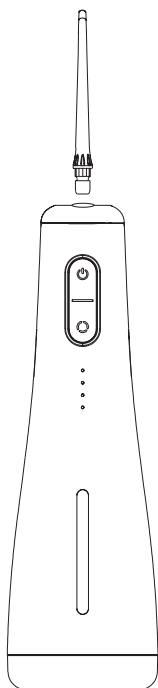KLAMER

Specifications

Model	KLAMER Water flosser
Mains voltage	3,7 V / 50-60 Hz
Power	8 Watt
Article no.	91476

Features

Thank you for purchasing the KLAMER water flosser. Please read this manual carefully BEFORE using the water flosser to ensure safe and effective use. The KLAMER water flosser is your ideal partner for daily oral hygiene and tooth cleaning.

- The water flosser provides effective and thorough cleaning of your teeth and gums.
- With four different modes (Clean, Soft, Pulse and DIY), you can customise the intensity of cleaning to suit your individual needs.
- Two different nozzles are included in the scope of delivery to provide you with the best possible support for your daily tooth cleaning.

Important Safeguards

Warning

The following basic safety precautions should be observed when using the KLAMER water flosser:


1. This appliance is not intended for use by persons (including children under 8 years of age) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they are supervised.
2. Cleaning and user maintenance must not be carried out by children without supervision.
3. Check that the voltage specified on the device matches the local mains voltage before connecting the device.
4. Do not use the device if the charging cable or the device itself is damaged. Always ensure that the plug is correctly inserted into a socket.
5. To protect against electric shock, the device must not be immersed in water or other liquids. Before each use, make sure that the cover for the charging port on the back is tightly closed.

6. Do not use any accessories that are not recommended by the manufacturer.
7. Do not modify or repair the product yourself if it is damaged and do not use it. Do not disassemble the product unless it is being disposed of.
8. If you have just undergone oral or dental surgery, are undergoing treatment or have other restrictions (such as weak teeth, fillings, bleeding gums, etc.), consult your dentist before using this product.
9. This product is intended for cleaning the mouth only and should not be used for any other purpose. If you experience discomfort or pain during or after use, stop using the product and consult your dentist.
10. If you continue to experience bleeding gums or other complications after two weeks of regular use, discontinue use of the product and consult your dentist.
11. Do not fill the tank with water at a temperature above 40°C.
12. Do not use this product without a nozzle fitted.
13. Please clean and dry the product after each use, otherwise bacteria and other pathogens may form on the surface.
14. This product contains a built-in rechargeable battery. Do not expose the device to particularly high or low temperatures and do not charge it in a damp or wet environment.
15. This product should only be used indoors, must not be exposed to direct sunlight and should be used at room temperature, otherwise malfunctions and failures may occur.
16. If battery fluid leaks and gets into your eyes, rinse them immediately with tap water and avoid rubbing your eyes. In case of complications, please contact your family doctor.
17. If battery fluid leaks and comes into contact with your skin or clothing, wash the affected areas thoroughly with tap water.

Product diagram

Explanations


- | | | | |
|---|---------------------------------|----|----------------------------------|
| 1 | Standard nozzle | 6 | Mode indicator/Battery indicator |
| 2 | Button for releasing the nozzle | 7 | Charging connector |
| 3 | Nozzle connection | 8 | Water tank filler opening |
| 4 | On/Off button | 9 | Water tank opening button |
| 5 | Mode button | 10 | Water tank lid |

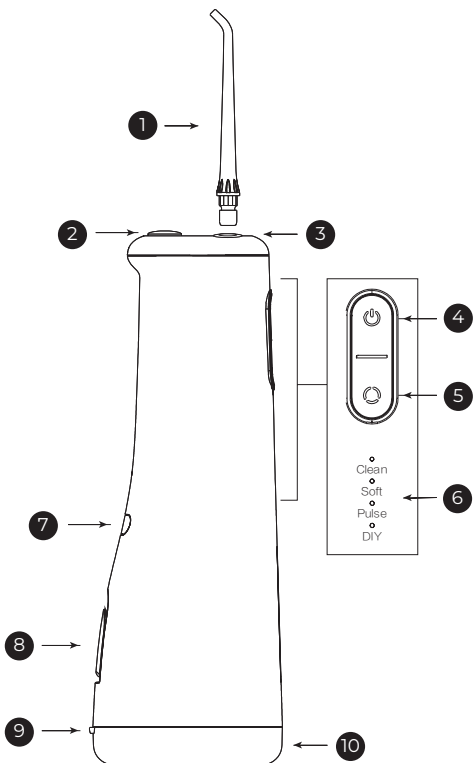Product diagram

Accessories

- 1 Charging cable
- 2 "Sensitive" nozzle

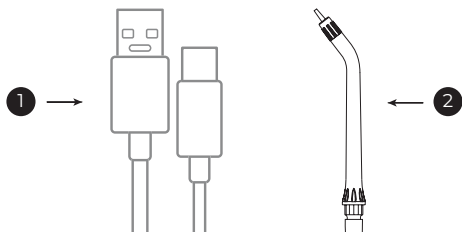
Note

If you have any questions or queries, please contact the KLAMER customer support team at hello@klamer-home.com. We will be happy to help you!

Tips and replacement parts

You can attach the second, unconnected nozzle to the travel holder in the water tank (❶) so that it is safely stowed away. The travel holder is located on the top of the water tank (❷).

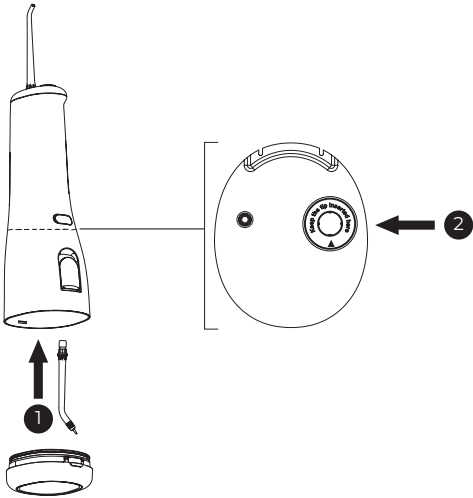Important

When closing the water tank with the lid, always ensure that the tank cap button is on the back of the appliance. Otherwise, the tank will leak and water may escape.

If the nozzles are lost or worn off, you can order replacement nozzles via the following link. Simply visit the website below or scan the QR code next to it.

Munddusche.klamer-home.com

Before first use

1. Take the KLAMER water flosser out of the packaging and remove all packaging materials.
2. Clean the appliance and the supplied nozzles with a damp cloth or, if necessary, with a little soap.
3. Fully charge the water flosser with the enclosed charging cable until all four points of the battery indicator light up.
4. Fill the water tank with clean tap water.
5. Place one of the supplied nozzles on the water flosser. Select the nozzle according to your individual needs (see p. 25) or as recommended by your dentist. Check that the nozzle is firmly inserted.
6. Carry out a test rinse to rinse out any production residues.

Operating the water flosser


To ensure the optimum performance and longevity of your KLAMER water flosser, it is essential to use and maintain it correctly. Please follow the steps below to ensure proper use:

1. **Charging the battery:** Before use, check that the battery is sufficiently charged. A new device must be charged for at least 4 hours before use. During the charging process and after switching off, the four dots on the battery indicator light up and show the battery charge status.


four lights = fully charged

one light = insufficiently charged

- 2. Nozzles:** Insert one of the two nozzles into the nozzle connection (❶) until you hear a click sound. To replace the nozzle, first switch off the appliance (❷), then press the button to release the nozzle (❸) and pull out the nozzle (❹).

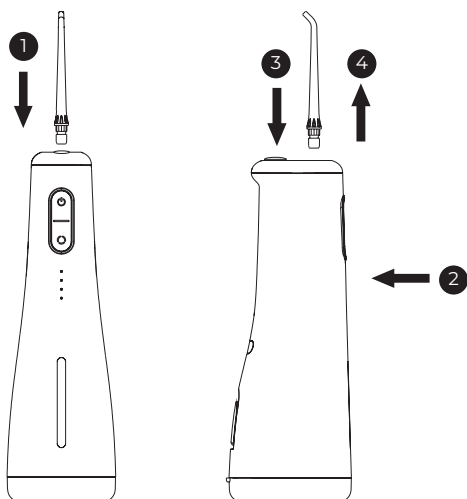
- 3. Water tank:** Fill the water tank with clean tap water or, if preferred, with mouthwash-enriched water. To do this, open the water filling connection, hold the water flosser horizontally (❺), fill the water tank (max. 250 ml) and then close the cap.

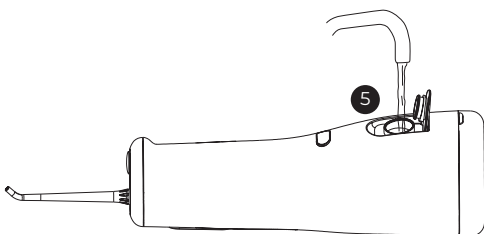- 4. Switch on the water flosser:** Press the On/Off button to switch the device on. Then select the desired mode using the mode button (an overview of the different modes can be found on (p. 25).

5. **Correct use:** Make sure that the nozzle is perpendicular to the cleaning area (6) and that the mouth is slightly open (7).

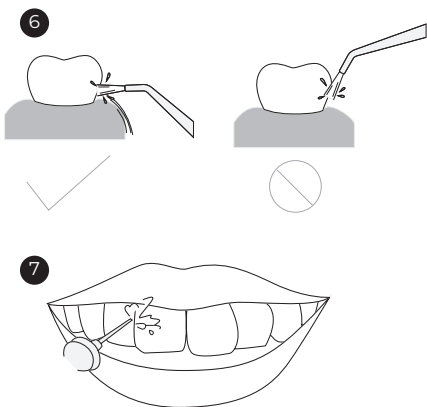6. **Cleaning process:** Start cleaning the desired areas. You can also use the Mode button during the cleaning process to change the cleaning mode to your own preference.
7. **After use:** After use, press the On/Off button to switch off the appliance. Pull out the nozzle as described in point 2, rinse the nozzle and empty the water tank completely. Dry the nozzle and the water flosser with a cloth.

The different modes

The KLAMER water flosser has four different modes, which differ in frequency and water pressure, so that you can customise the cleaning to your individual needs.

Mode	Water pressure	Explanation
Clean	ca. 860 kPa	Standard mode for daily cleaning.
Soft	ca. 590 kPa	Lower water pressure for sensitive gums.
Pulse	approx. 460 kPa and 860 kPa	Mode with alternating water pressure.
DIY	from ca. 340 kPa to 860 kPa	Press and hold the mode button to set the preferred water pressure.

The different nozzles

The scope of delivery includes two different nozzles for the KLAMER water flosser, which provide you with the best possible support for your daily cleaning.

Nozzle	Explanation
Standard nozzle	Effectively removes food debris from the teeth, gums and interdental spaces.
Sensitive nozzle	Has a softer, smaller mouthpiece than the standard nozzle, which allows for gentler cleaning.

Note

When using for the first time or with sensitive gums, we recommend using the “Soft mode” in conjunction with the “Sensitive nozzle” in order to gently introduce the gums to the cleaning process.

Cleaning

Before cleaning, make sure that the charging connection is firmly closed, the electronics in the motor element should NEVER come into contact with water. We recommend that you clean the water flosser completely at least once a week.

You clean the KLAMER water flosser as follows:

1. Remove the water tank cover and empty the water tank. Use a slightly damp cloth to clean the tank. Then dry it thoroughly.
2. Press the on/off button so that the water flosser pumps the excess water out of the hose and the nozzle. Then make sure that the appliance is switched off.
3. Clean the outside of the water flosser with a slightly damp cloth. Make sure that no moisture penetrates the appliance.
4. Do not use any aggressive cleaning agents or abrasive cleaners for the water flosser to avoid damaging the surface and the sensitive nozzles.
5. Remove the attached nozzle and rinse both nozzles under running water to remove any residue. For more thorough cleaning, the nozzles can be immersed in a mild vinegar solution (equal parts water and vinegar) to dissolve limescale deposits.
6. Allow all parts to dry completely before reassembling them and using the water flosser again.

Note

If you have used the water flosser with a mouthwash-enriched mixture, residues may form in the tank. These areas must be cleaned more thoroughly.

Storage

Make sure that the water flosser, especially the water tank and the nozzles, are clean and dry. Follow the cleaning instructions in the “Cleaning” section.

Store the device as follows:

1. Store the water flosser in a cool, dry place out of direct sunlight. To protect the water flosser as much as possible, we recommend storing it in its original packaging or another suitable container that protects it from dust and moisture.
2. Avoid placing heavy objects on the water flosser or accessories to prevent damage to the sensitive parts.
3. Make sure that the appliance and its components are completely dry before storing them. This helps to prevent the formation of mould and bacteria as well as possible malfunctions due to moisture.
4. As already mentioned on page 21, you can place the second, unconnected nozzle on the travel holder in the water tank so that it is safely stowed away.

Problem solving

In the following table you will find a list of frequently occurring problems, their causes and solutions that can help you during use:

Problem	Cause	Solution
The device cannot be switched on.	Device purchased new or not used for more than three months.	Charge the device for 4 hours so that all 4 points of the charging indicator light up.
	The device is out of order.	Please contact the KLAMER customer support team at hello@klamer-home.com .

Note

Please contact the KLAMER customer support team at **hello@klamer-home.com** if you have any further questions or queries. We will be happy to help you!

The water pressure is too low during use.	A foreign object in the nozzle is blocking the water outflow.	Clean the nozzle.
	A foreign object in the filter tip of the water tank hose is blocking the water flow.	Clean the filter tip.
The battery discharges after a very short time.	The battery has reached the end of its useful life.	Please contact the KLAMER customer support team at hello@klamer-home.com .
	The battery has not been charged long enough.	Charge the device for 4 hours so that all 4 points of the charging indicator light up.
The charging indicator does not light up when charging.	The charging cable or the mains plug are not connected correctly.	Make sure that both the charging cable and the mains adapter are plugged in correctly.
	The charging circuit of the product is defective	Please contact the KLAMER customer support team at hello@klamer-home.com .
	The original charging cable is not used.	Use the original charging cable.
Water is leaking from the water tank.	The lid was not placed on the water tank properly.	Please follow the instructions on p. 21.
	The lid was damaged.	Please contact the KLAMER customer support team at hello@klamer-home.com or order spare parts from the following website: munddusche.klamer-home.com .

Warranty

The manufacturer provides a warranty in accordance with the legislation of the customer's country of residence for a minimum period of 12 months, starting from the date on which the device is sold to the end user.

The warranty covers only defects in materials or workmanship, which does not apply in cases of defects in materials or workmanship:

- Normal wear and tear
- Improper use, e.g. overloading of the device, use of non-approved accessories
- Use of force, damage due to external influences
- Damage due to non-observance of the user manual, e.g. connection to an unsuitable power supply or non-observance of the cleaning instructions
- Partially or completely dismantled devices

Conformity

This product conforms to the applicable CE directives. The full Declaration of Conformity can be requested using the contact form on our website or by contacting KLAMER Customer Service at the e-mail addresses given in this manual.

Notes on disposal

If the image below (crossed-out wheeled garbage can) is on the product, European Directive 2002/96/EC applies. These products must not be disposed of with normal household waste. Inform yourself about the local regulations for separate collection of electrical and electronic equipment. Follow the local regulations and do not dispose of old devices with household waste. Disposing of old equipment in accordance with the regulations will protect the environment and the health of your fellow human beings from possible negative consequences. Material recycling helps to reduce the consumption of raw materials.

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