Warranty

The manufacturer provides warranty in accordance with the legislation of the customer's own country of residence, with a minimum of 12 months, starting from the date on which the appliance is sold to the end user.

The warranty only covers defects in material or workmanship, which will not apply in cases of:

- Normal wear and tear
- Incorrect use, e.g. overloading of the appliance, use of non-approved accessories
- Use of force, damage caused by external influences
- Damage caused by non-observance of the user manual, e.g. connection to an unsuitable mains supply or non-compliance with the installation instructions
- Partially or completely dismantled appliances.

Declaration of Conformity

This product is in accordance with the following European directives:

2014/30/EU (EMV)

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Disposal Information

If the icon below (crossed-out trashcan on wheels) is found on the product, the European directive 2002/96/EG applies. These products may not be disposed of in the household waste at the end of its useful life. Disposal can take place at appropriate collection points provided in your country.

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
DEEP TISSUE NECK MASSAGER

USER MANUAL

Read this manual thoroughly before using and save it for future reference

Fitleap
Deep Tissue Massager

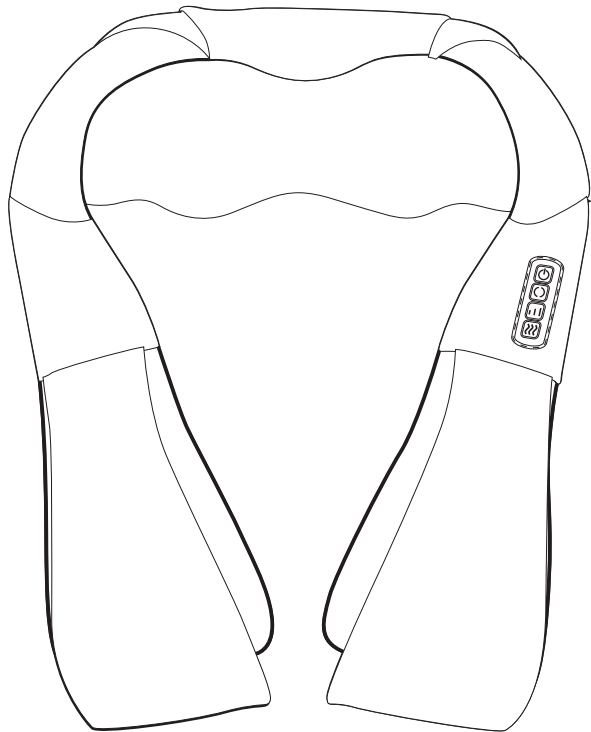Find help and advice at
neckmassage.fitleap.de
or scan the QR code on the left

Technical Details

Model	Fitleap Deep Tissue Neck Massager
Rated Voltage	220-240 V / 50-60 Hz
Rated Power	30 W
Article No.	63859

Note

Please contact the Fitleap customer support team if you have any questions or problems. You can reach us via hello@fitleap.de. We're happy to help!

Features

Thank you for purchasing the Fitleap Deep Tissue Neck Massager. For your own safety, please read this manual carefully BEFORE using the device. There are many benefits to using the neck massager:

- The neck massager is especially helpful for tension and pain in the neck and shoulder area.
- The neck massager can be used to massage lactate out of muscle tissue in case of sore muscles
- The neck massager can also be used before or after sports to relieve tension
- Different heat and massage intensities help to find exactly the right combination for the muscle

Important safety precautions

Warning

Before using the neck massager on your body, please get your doctor's approval first if any of the following factors apply to you:


- Pregnancy, diabetes with complications (e.g., neuropathy or retinitis), wearing a pacemaker
- Peripheral neuropathy, recent surgery, epilepsy, migraines
- Herniated disc, spondylolisthesis, spondylolysis, or spondylosis,
- Recent joint replacement or IUDs, metal pins or plates, or other concerns about your physical health and constitution, especially in the neck, shoulders, or throat.

In general, the fitness gun is not suitable for use by children and persons under 13 years of age.

Important

1. This device is not intended for use by persons (including children under 13 years of age) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they are supervised. Keep the device out of the reach of children.
2. Do not use the device as a toy. Special attention is required when using near children.
3. Never immerse the device in water or other liquids. Never pour water over the device.
4. Avoid sources of heat and direct sunlight to prevent the device from heating up excessively.
5. Do not use the device for any purpose other than personal massage.
6. Do not touch any part of the power supply with wet hands. Never unplug the power supply by pulling on the cord, always pull on the plug or power supply itself.
7. Use the device only with the supplied power supply.
8. Do not use lubricants, cleaners, polishes or air fresheners on this equipment. Contact Fitleap support team if repair is required. Never disassemble the unit by yourself, as improper assembly may result in electric shock or fire.
9. If the device is defective, has been dropped from a great height, damaged, left outdoors, received a violent shock or fallen into water, it should only be put back into operation carefully and the Fitleap support team should be consulted if necessary.
10. The device works best at normal room temperature
11. Use the neck massager while driving only as a passenger, not as a vehicle driver.
12. Do not use this product if you have applied oils or lotions to your skin, or if you have sweated profusely.
13. Do not operate the device under a blanket or pillow. Excessive heating may cause fire, electric shock or injury.
14. We recommend not to use the device before going to sleep. The massage may have a stimulating effect and delay sleep.

Product diagram

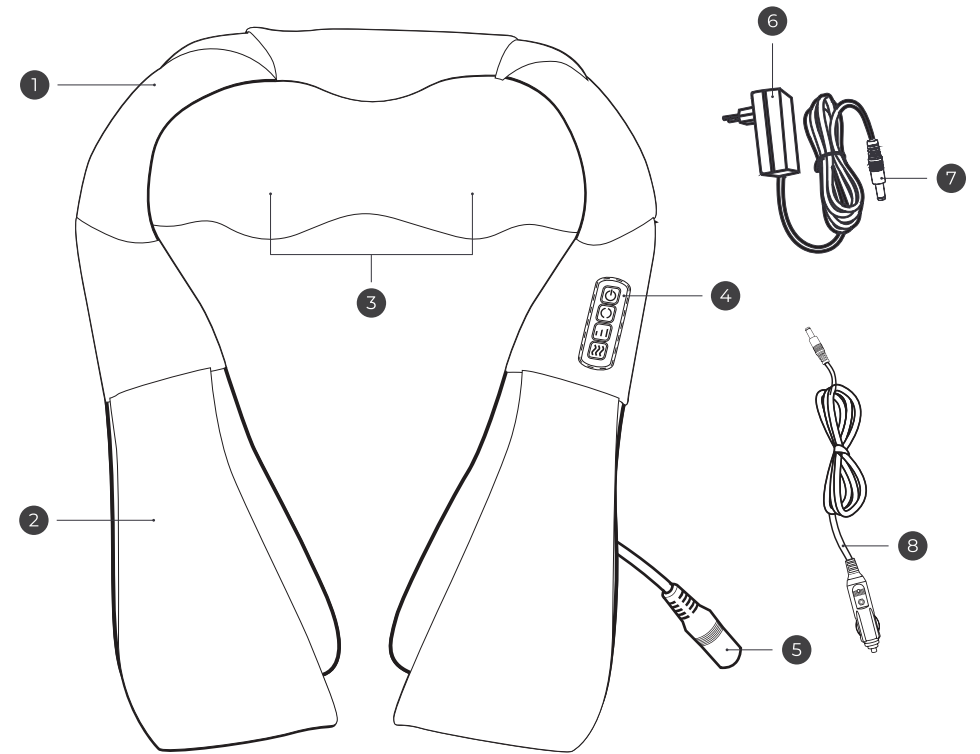- 1 Leather sheath
- 2 Arm loops
- 3 Massage heads
- 4 Control element

- 5 Power cable connection
- 6 Power supply
- 7 Power cable
- 8 Car power supply with cable

About neck massage and shiatsu

The Fitleap neck massager is based on the principles of the traditional Japanese massage technique Shiatsu.

The term comes from Japanese and means "finger pressure". It is not a massage in the classical sense. The basis of Shiatsu massage is the meridian system of the body, i.e. the energy channels through which life energy flows.

Pressure, touch and manipulation techniques are used to take a natural approach to promoting mental and physical well-being.

Shiatsu massage can help with insomnia, digestive problems, poor posture, headaches, coughs and colds, sprains and especially tension in the neck, shoulders and back.

The buttons of the control element

The neck massager is operated using the four buttons on the control element. This is located on the outer side of the left arm strap and have the following functions:

Symbol	Name	Meaning
	Power button	Pressing the power button turns the device on and off.
	Speed button	The speed button can be used to select three different speed levels with increasing intensity.
	Direction button	The direction of the massage can be selected via this button. It changes automatically every minute.
	Heating button	The integrated heat function of the massage element can be controlled via the heat button. It is automatically active when the massage is started or the device is switched on.

Using the neck massager

1. Connect the device to the power supply by plugging the power supply into a suitable socket (or in the case of mobile use, into the cigarette lighter of your car).
2. Position the device in the desired location or position. For the neck massage, put your arms through the hand loops and apply a little counterpressure to ensure that the device remains in position on the neck.
3. Now switch on the device by pressing the power button on the control element. The device will start massaging and the area behind the massage heads will light up red.

Setting the massage speed

The massage speed has three different levels. When starting the massage, the lowest level is selected. By pressing the button you switch to the next higher level. If the highest level is already selected, then the device switches back to the lowest level.

Changing the massage direction

The massage direction changes by itself every minute. If you want to change the direction manually, then this is possible by simply pressing the corresponding button. Pressing it again will change the massage direction back to the original direction.

Turning the heating function on / off

The heat function of the neck massager is automatically activated when the massage is started or the device is switched on. It can be deactivated by simply pressing the heat button and reactivated by pressing it again.

Using the massager for different muscle groups

You can generally use the neck massager on the whole body. Possible restrictions and notes on non-recommended use can be found on page 10 under *Important safety precautions*.


Below you will find some examples for the application on different massage zones or muscle groups. You can adjust the pressure of the massage using the arm loops.

Examples of different muscle group applications


Be careful not to apply too much pressure so as not to damage the device. If the massage heads slow down noticeably or stop completely, the pressure is too strong and you should apply less pressure.


Neck & Shoulder

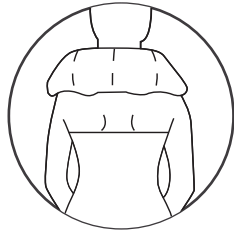
Neck

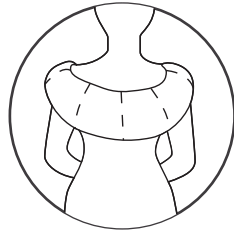Upper Back

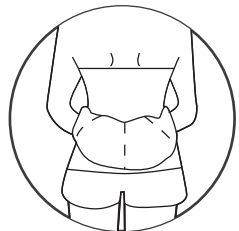
Lower Back

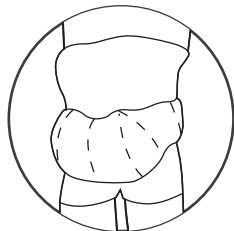
Abs


Hip

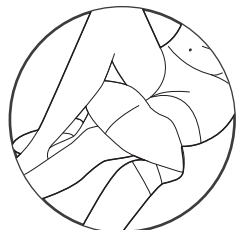Legs

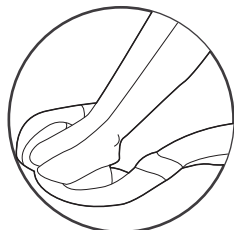Feet

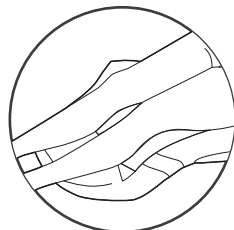Thighs

Note

The neck massager switches itself off after 15 minutes of continuous use to prevent overloading or overheating of the device. If you massage different massage zones one after the other, allow the device to cool down for a few minutes in the meantime. This applies in particular if you also use the heat function.

Cleaning and storage

Never clean the device while the power supply is connected to the socket and make sure that the device is switched off when cleaning it. This will prevent the risk of electric shock.

Do not use water or harsh or corrosive cleaners or solvents for cleaning. Instead, take a slightly damp cloth and gently wipe the individual components of the neck massager (except the power supply). Dry them carefully afterwards.

Problem solving

Problem	Possible cause	Solution
The neck massager can't be turned on.	The device does not have power.	Make sure that the power supply is connected to a power source and that the cables are connected correctly.
	The device is damaged or broken.	Contact the Fitleap customer service
The neck massager suddenly turns off during use.	The automatic timer of 15 minutes has activated to avoid overheating.	Allow the device to cool down for a few minutes before continuing the massage.
	The power supply is no longer connected correctly or the cables are no longer connected.	Check the power supply to the device.
	The device is damaged or broken.	Contact the Fitleap customer service
The warming effect of the heating function can't be felt / is very light.	You wear clothes that do not conduct heat well.	Wear fewer and less thick / dense layers of clothing
	The heating mode is not switched on.	Switch on the mode using the corresponding button.
	The heat mode or the heat mode button is defective	Contact the Fitleap customer service
The massage is too intense / not intense enough	You are applying too much / too little pressure for the current massage zone via the hand loops.	Control the pressure according to your needs by applying less or more tension to the arm loops.

Note

Please contact the Fitleap customer support team if you have any questions or problems. You can reach us via hello@fitleap.de. We're happy to help!